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NEU
NEW
250.000
Lichtimpulse
Light pulses

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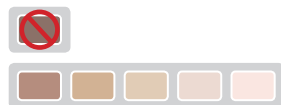
Art.-Nr.: 576.09

IPL 10000+ SalonPro System

Technology for long lasting hair removal used in professional cosmetic institutes**

Safe, fast and reliable

A lifetime of flashes - includes cartridge with 250,000 light pulses****

Fastest Application - Around just 26 min. to treat the entire body***

XXL 7cm² light surface for especially fast treatment

Automatic skin tone detection

Clinically tested technology for maximum safety at home

Maximum safe application with the 2-in-1 skin-type and skin-contact sensor

Incl. precision attachment

Up to 50% hair reduction after 3-4 treatments

State-of-the-art IPL light technology used by dermatologists for professional hair removal. Deactivates the hair roots deep beneath the skin and in doing so decreases the chance of the hair growing back.


With "flash-glide" function for comfortable usage

Integrated UV filter to protect the skin

5-year warranty

Cannot be used on dark skin

May be ineffective on white or light-coloured hair

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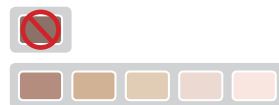* Dermatological study, 36 women, Institut proDERM Hamburg, study no. 14.0148-11
 ** Technology appraised by Institute proDERM, Hamburg
 *** On the highest power setting for arms, legs, face, armpits and the bikini line
 **** Sufficient for 40-60 years with normal usage

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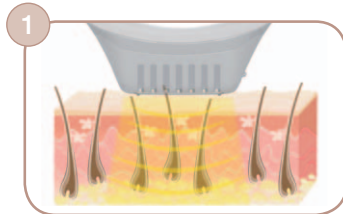
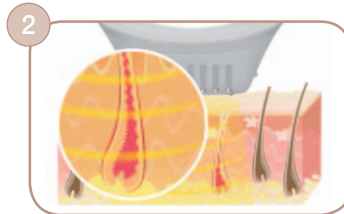
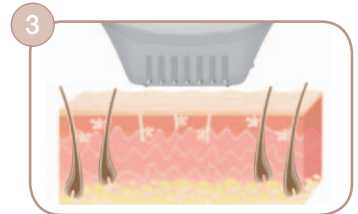
Skin type \ Hair	I	II	III	IV	V	VI
Black	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Dark brown	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Brown	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Dark blonde	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Light blonde	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Red head	Energy 4-5-6	Energy 4-5-6	Energy 3-4	Energy 1-2-3	Energy 1-2-3	✗
Gray/white	✗	✗	✗	✗	✗	✗

* Energy should primarily chose based on your own comfort
** The higher the energy, the more efficient the treatment will be

Anwendung | Application

Körper | Body

Gesicht | Face